


# ACTU8 TALON™

Dynamic Pressure Controlled Forceps

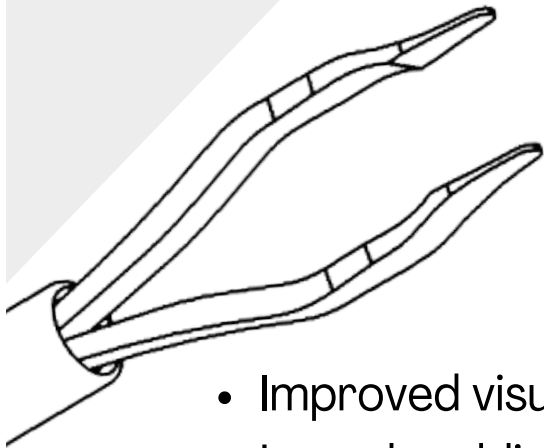


## Your Next Generation ACTU8 Forceps!

- Improved visualization
- Less shredding of membranes
- Increased grams of force
- All of the benefits of the ACTU8 handle
- Universal applications including ILM, ERM PVR, dropped IOLs & Diabetic Membranes

### What is DPC??

Dynamic Pressure Control gives the control of the forceps back to the surgeon. As the surgeon depresses the forceps, it is not an all or nothing with grip strength. Talon tips spread the pressure throughout the body of the tip. This allows for a more gradual build up of grip strength leading to less shredding of membranes and more surgeon control.

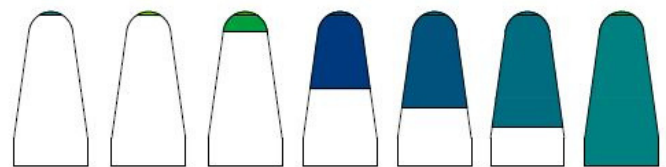

TALON TIP PRESSURE AS THE FORCEPS CLOSES

### Responsive.

- Linear Actuation
- Instantaneous Tip Movement
- Low Closure Force

### Consistent.

- Patented Tip Fabrication Process
- Robust Inspection Procedure
- Reproducible Results, Every Time

### Effective.

- Excellent Maneuverability
- Superior Grip Strength<sup>1</sup>
- Continuous Evolution

MKT-0001 rev A

Catalog Numbering and Ordering information are located on the reverse side

[www.vortexsurgical.com](http://www.vortexsurgical.com)

[info@vortexsurgical.com](mailto:info@vortexsurgical.com)





# Catalog Numbers and Descriptions

## TALON PRO-GRIP<sup>2</sup>

VS0751.23 23ga Talon Pro-Grip Forceps  
VS0751.25 25ga Talon Pro-Grip Forceps

## TALON ILM

VS0752.23 23ga Talon ILM Forceps  
VS0752.25 25ga Talon ILM Forceps

## TALON END-GRIP

VS0753.23 23ga Talon End-Grip Forceps  
VS0753.25 25ga Talon End-Grip Forceps

## ORDERING INFORMATION:




**ACTU8TALON**  
Dynamic Pressure Controlled Forceps

Vortex Surgical  
4 Research Park Drive, Suite 124  
Saint Charles, MO 63304 USA

Phone: 636-778-4350  
Email: [info@vortexsurgical.com](mailto:info@vortexsurgical.com)  
Web: [www.vortexsurgical.com](http://www.vortexsurgical.com)

1-Vortex Data on file 2-Coming Soon

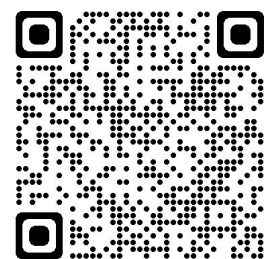

**VORTEX SURGICAL**